
CREATE

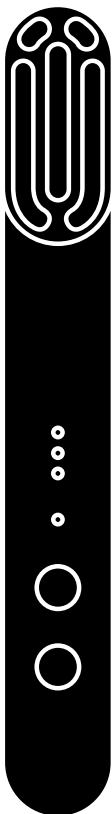Skin Care Led Wand 4 in 1

User manual | Manual de instrucciones

[illegible]

[illegible]

CREATE

INDEX

ENGLISH

Safety instructions	6
Product Features	7
Parts list	8
Button/Mode Description	8
Instructions for use	9
Battery	10
Technical specifications	10

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ENGLISH

Thank you for choosing our product. Before using this device and to ensure its best use, please read the instructions carefully.

This manual contains safety information for the correct use of the product: read it carefully and retain it to avoid damage and hazards caused by improper use. The product is for personal use and should not be used by multiple people.

SAFETY INSTRUCTIONS

Before using this product, always read the basic safety precautions.

- Do not disassemble or modify this product without authorization.
- Do not turn the head in the opposite direction to that indicated in the instructions.
- This device is not waterproof: do not wash or submerge it.
- When used near the eyes, a flickering sensation may be felt; this is a normal bioelectrical effect of microcurrent.
- Since the body's electrical resistance varies from person to person, adjust the intensity according to your tolerance and discontinue use immediately if discomfort occurs.
- After each use, clean with a dry cloth or slightly damp wipe.
- If you notice any abnormalities (such as no vibration), check the battery level or if the device has turned off.
- When not in use for a long period, store it in a dry place away from direct sunlight.
- Keep this product out of reach of children; it is not a toy.
- **WARNING!** This product should not be used by people with electronic medical devices such as pacemakers or other implants; with heart disease, malignant tumors, or irregular blood pressure, with infectious diseases or acute illnesses accompanied by fever, during dental treatment, with weakened immunity or under immunosuppressive therapy, with hemophilia, coagulation disorders, diabetes, facial neuralgia, or changes in sensitivity and pain, with varicose veins, facial melanoderma, light allergy, or abnormal pigmentation, with microvascular dilation caused by prolonged use of steroid hormones or liver dysfunction, with skin allergies, atopic dermatitis, skin diseases,

or very sensitive skin, with a history of epilepsy or discomfort caused by flashing lights, with after-effects on the head or eyes, or who have undergone eye surgery (in this case, consult your doctor before use), pregnant or breastfeeding women; young children, people on medication that causes drowsiness, with extreme fatigue, or in poor physical condition, people with skin problems related to pigmentation, inflammation caused by chemical fibers, or adverse reactions to cosmetics, as well as those who have recently undergone chemical peels or other cosmetic procedures that have caused irritation or injury.

- Do not use on parts other than those designated.
- Patients under medical treatment should consult their physician before using the product.

PRODUCT FEATURES

EMS MICROCURRENT

Low-level electrical current that stimulates the face by sending gentle waves through the skin to the facial muscles. It tones muscles, stimulates their vitality and promotes the production of collagen and elastin.

HEAT FUNCTION

This function features a constant 42°C intelligent temperature control system that relieves facial fatigue, promotes blood circulation, and enhances the absorption of creams and cosmetics. Continued use can help resolve facial problems such as dark circles, bags, and puffiness.

OPTICAL BEAUTY

Its 630 nm red light is absorbed by the dermis, accelerating cell regeneration and promoting collagen synthesis, thereby achieving a filling effect on epidermal cells, firming the skin, fading acne marks and pigmentation, and eliminating fine lines.

SONIC VIBRATION

When used with eye serum, eye cream, or facial care products, sonic vibration helps break down the essence into tiny molecules, promoting its absorption.

MEMORY MODE AND INTELLIGENT TIME CONTROL


Each time it is turned on, it memorizes the last mode used, saving time and making it easier to use. Additionally, its automatic time control feature turns off the device after 5 minutes of use to prevent unintentional prolonged use that could damage the skin.

90° ROTATING HEAD

The innovative 90° rotation allows the product to better adapt to different areas, such as the eye and lip areas, as well as the skin on the entire face and neck.

PARTS LIST

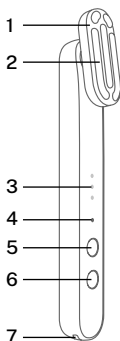1. Light emission zone
2. Head
3. Intensity/charge indicator light
4. Mode indicator light
5. EMS intensity adjustment button
6. Power/Mode Change Button
7. Loading port

BUTTON/MODE DESCRIPTION

EMS INTENSITY ADJUSTMENT BUTTON

Three skincare modes, each with three intensity levels. The intensity setting is suitable for different skin types depending on their tolerance level:

- The first level is suitable for sensitive skin.
- The second level is suitable for normal skin.
- The third level is suitable for more resistant skin.

The three intensity indicators are white lights, corresponding to the three levels.

POWER BUTTON / MODE SWITCH


Long press the power button to turn it on, short press to adjust the level and adjust the mode to get different effects:

- **Relaxing mode:** soothing microcurrent + blue light + vibration. Promotes blood circulation and relaxes facial muscles. The second level is suitable for normal skin.
- **Firming mode:** contraction microcurrent + red light + vibration. Effectively contracts muscles, with a lifting, firming, and swelling-reducing effect.
- **Activation mode:** microcurrent activation + purple light + vibration. Stimulates ATP production and promotes collagen regeneration.

INSTRUCTIONS FOR USE

It is recommended to use it after facial cleansing to prevent facial oil from affecting the user experience.


1. Pay attention to the direction of the rotation arrow on the product. Simply turn the head clockwise from the “I” state to the “T” state, and then turn it counterclockwise again. Do not rotate the head in the opposite direction, as this will cause damage to the head by external force and affect its operation.

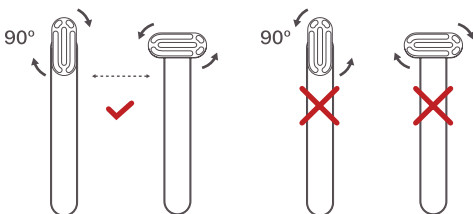2. When using the device near your eyes, light flickering may occur due to the microcurrent function and is a completely normal phenomenon; you should not worry about it. For best results, apply eye contour or facial cream evenly. If you're using it in winter, you can start by applying it to your hands for a few seconds, wait for it to warm up, and then use it on your face.

Note: When the battery runs out of power, the device beeps twice and switches off automatically.

3. For facial use, continuous application can effectively solve skin problems such as fine lines around the eyes, nasolabial lines, neck lines, etc., as well as firm and smooth the skin. For different areas of the face, the “I” and “T” states can be adjusted to appropriate angles according to each area. The method of use is shown in the figure below. It is recommended to apply it from the bottom up to achieve the lifting effect.

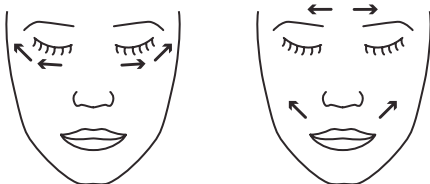BATTERY

Battery information:

- When the battery is low, the three white lights flash rapidly.
- During charging, the three white lights turn on alternately.
- When charging is complete, the three white lights remain on.
- Charging time: 120 minutes. Duration: 60 minutes.
- When the battery is exhausted, please charge it in time.

LOADING METHOD

Insert one end of the USB charging cable into the bottom charging port of the product and connect the other end to the charger or a computer's USB port. You will know the product is charging correctly when the three white lights flash alternately.

The recommended time for each charge should not exceed 4 hours.

TECHNICAL SPECIFICATIONS

PARAMETER	WORTH
Power	1.5 W
Strain	3.7 V
Battery	420 mAh
Frequency	Mode 1: 10 Hz Mode 2: (5 Hz-20 Hz-5 Hz)/9 s Mode 3: (20 Hz-333 Hz-20 Hz)/6 s
Cable type	USB-C
Loading time	<120 min
Operating time	>60 min
Automatic shutdown	In 5 minutes
Dimensions	160.5 x 21.8 x 25.7 mm
Weight	65 g
Materials	Stainless Steel + Aluminum Alloy + PC + ABS

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable re-use of material resources. To return your used device, please use the return and collection systems or contact the retailer from where the product was purchased. They can take this product for environmental safe recycling.

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