CREATE

Silk Skin IPL

User manual | Manual de instrucciones

INDEX

ENGLISH

Safety instructions	6
Product Introduction	8
Parts list	8
Skin sensitivity test	8
How to use	9
Course of treatment	11
Replacing the lamp	11
Skin rejuvenation function	11
Areas suitable for treatment	12
Adverse reactions and related countermeasures	12
Technical specifications	12
Troubleshooting	13

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ENGLISH

Thank you for choosing our pulsed light epilator. Before using this appliance and to ensure its best use, please read the instructions carefully.

The safety measures outlined here reduce the risk of fire, electric shock, and injury when followed correctly. Keep the manual in a safe place for future reference, as well as the warranty, purchase receipt, and box. If applicable, give these instructions to the future owner of the appliance. Always follow basic safety instructions and risk prevention measures when using an electrical appliance. The manufacturer shall not be liable for any damage resulting from the user's failure to follow these instructions.

SAFETY INSTRUCTIONS

When using any electrical appliance, the following basic safety precautions should always be observed.

- The product may only be used for its intended purpose. The manufacturer shall not be liable for damage resulting from improper use or incorrect handling.
- This device contains high voltage. Do not disassemble the device.
- Do not block the vents during use. They are used to cool the product.
- Do not allow water or moisture to enter the equipment.
- Do not use the device if the plug or cord is damaged.
- The use of the device in flammable or explosive conditions is prohibited.
- Do not use the product under any of the following conditions:
 - In children and the elderly.
 - People with local or systemic inflammation; people with photosensitive skin or those using photosensitive medications.
 - Users of implanted electronic medical devices (e.g., cardiac pacemakers) and life support systems (e.g., artificial heart-lung machines), as these may cause malfunctions of the electronic medical devices.
- The product is not recommended for use by women who are pregnant, attempting to become pregnant, or breastfeeding, as this product has not been tested in these groups.
- People who have sunbathed, used tanning beds, or used artificial tanning products should not use the product immediately. If you are unsure whether your skin has been tanned or sunburned, it is recommended to perform a skin sensitivity test and wait 24 hours.
- Do not use directly on tattoos or dark or black spots, as they can absorb large amounts of light and cause skin discomfort.

- Use on sensitive areas such as the eyes, neck, nipples, genitals, or anus is prohibited, as the skin in these areas is especially sensitive and vulnerable.
- Do not apply for too long to the same area or repeat the process too many times, as this may cause excessive heat and discomfort to the skin.
- If your skin appears red and swollen for more than 24 hours after use, discontinue use immediately and consult your doctor or customer service.
- Do not direct the lamp towards vulnerable parts such as the eyes when emitting light. Always wear eye protection during use.
- Do not look directly at the light emitted by the lamp or even at the reflected light, even if you are wearing eye protection.
- Before use, remove all metal objects such as rings or bracelets to prevent equipment malfunction.
- Keep the room ventilated.
- Keep the product away from children and store it out of their reach after use.
- It is strictly prohibited to disassemble or modify this product and operate it improperly. If disassembled without authorization, the product will immediately lose its right to warranty service.
- Apply aloe vera gel or a moisturizer to the treated skin. Do not use irritating cosmetics immediately.
- Make sure the treated skin is clean before starting the procedure.
- For 2 hours after treatment, do not wet the area with water. Do not sunbathe, take hot baths, steam baths, etc., for the next 48 hours.
- Avoid exposing treated skin to direct sunlight for 2 hours after treatment. It is recommended to perform the treatment at night and, if you go out, apply sunscreen.
- Do not swim in pools for 48 hours after treatment.
- If you ignore these warnings, using the device may cause skin redness and swelling, as well as temporary or permanent color changes.
- Do not press the light trigger button continuously. Just press it once to emit a flash. It is recommended to irradiate each area 2 to 3 times (with the lamp in full contact with the skin).
- Do not use the appliance continuously for more than 25 minutes. If the treatment area is large, divide it into several sessions. After 25 minutes of use, unplug the appliance, let it rest for 10 minutes, and then continue.

PRODUCT INTRODUCTION


This beauty care device is an innovative, multifunctional tool that combines laser hair removal, skin rejuvenation, acne treatment, and other benefits in a single device. Thanks to the integration of advanced optical technology, it presents itself as a professional, practical, and safe solution in the field of beauty and personal care.

Before using the device, it is essential to carefully read all instructions and precautions to ensure correct and safe use. Its design prioritizes both user comfort and protection, and its manufacturing meets the highest quality standards.

It is indicated for men and women who wish to effectively remove unwanted hair. The treatment is a completely physical and non-invasive therapy. For darker skin, it is recommended to start with the lowest intensity and gradually adjust the energy level according to skin tolerance and response.

PARTS LIST

1. Light output
2. Skin detection sensor
3. Light trigger button
4. Power port
5. Level indicator
6. On/off/level adjustment button
7. Rejuvenation lamp: Red filter
8. Hair removal lamp: Yellow filter
9. Ready indicator light

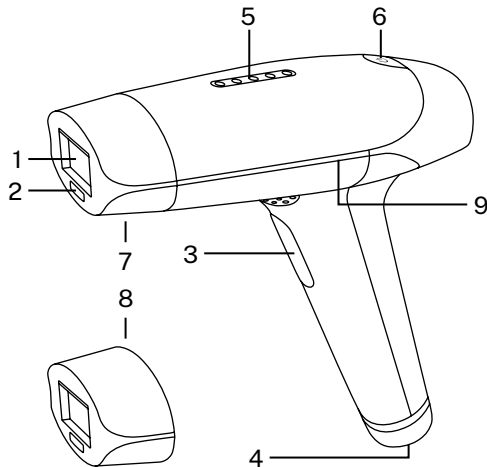SKIN SENSITIVITY TEST


If this is your first time using the product, or your first time applying it to a new area, first select a small area to perform the sensitivity test by following these instructions:

1. Turn on the device following the instructions shown below.
2. Choose a low level: level 1 or 2.
3. Perform the test in the selected area following the operating instructions and paying attention.
4. Check the condition of the skin after 30 minutes.
5. If the skin in the test area is normal, you can continue using the product. If the skin becomes red or swollen, wait 24 hours and repeat the test.

HOW TO USE

1. Before use, shave all hair around the treatment area. Clean and dry the area thoroughly to maximise the amount of energy absorbed by the hair follicles for the best treatment effect.

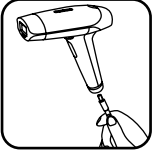2. Insert the plug of the power adaptor into the power socket of the appliance.

3. Plug the adapter into the power socket, the level indicator will light up for 2 seconds, this indicates that the device has been connected to the power correctly.

4. Press and hold the power button for 2 seconds to switch on the device.

5. The device has 5 levels (1-5). Level 1 is the default when you turn it on. Briefly pressing the level adjustment button increases the level by one level at a time, with the corresponding indicator lighting up. Once level 5 is reached, pressing it again returns to level 1, operating in a loop. Always start at a low level and gradually increase it as appropriate. Adjust the level according to your own skin's adaptability.

6. Use the glasses exclusively for this product to avoid eye damage.

7. Manual and automatic light emission settings

- a. **Manual light emission settings.**

The default setting is manual light emission, meaning the device will flash once each time the light trigger button is pressed (suitable for hair removal on small areas).

- b. **Automatic light emission settings**


After turning on the device, press and hold the light trigger button for 2 seconds. The readiness indicator light will flash and auto-broadcast mode will activate. You just need to bring the lamp close to the skin and the device will automatically emit light (suitable for hair removal on large areas).

- c. **Attention**

When the device is in automatic emission mode, when pressing the light trigger button you must be careful not to expose other areas to the lamp and never allow other people to touch it.

- d. **Switch to manual light emission**

When the device is in automatic emission mode, press and hold the light trigger button for 2 seconds. The readiness indicator light will stop flashing and the device will return to manual broadcast mode.

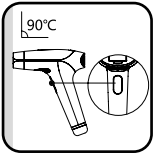8. Using the lamp in the hair removal area

Place the device in direct contact with the area to be treated and try to keep it at a right angle to the skin.

a. Manual emission mode:

When the device is in manual emission mode and the readiness indicator is flashing, press the light trigger button to irradiate. The area covered by the lamp will be hair removed. Then, move it to the next treatment area and, when the readiness indicator flashes again, press the light trigger button again to irradiate. Repeat these steps. It is recommended to irradiate each area 2 to 3 times.

b. Automatic emission mode:

When the device is in auto-emitting mode, the device will automatically emit light, so there is no need to press the light trigger button. The area covered by the lamp will be hair removed. Then, move it to the next treatment area and make sure it is in contact with your skin; the device will automatically irradiate. Repeat these steps. It is recommended to irradiate each area 2 to 3 times.

Advice: When the device is in automatic emission mode and the lamp remains fixed on the skin without moving, it will irradiate every 8 seconds.

9. After use, press and hold the power button for 2 seconds to turn off the device. It will automatically shut off in 5 minutes if no action is taken.
10. Unplug the adapter and remove it from the power outlet.
11. After 10 minutes, when the lamp has cooled, clean the lamp's light outlet with a damp, alcohol-free tissue or soft tissue.
12. Once clean, store the device and adapter in the package box. Place the box in a dry, well-ventilated place, out of the reach of children.
13. Care of treated skin: After 20 minutes, clean the treated skin. You can apply moisturizing or repairing products to nourish your skin. If you experience discomfort, apply cold to the area for 30–60 minutes until the burning sensation disappears.

ATTENTION

- Apply aloe vera gel or a moisturizer to the treated skin. Do not use irritating cosmetics immediately.
- For 2 hours after treatment, do not wet the area with water. Do not sunbathe, take hot baths, steam baths, etc., for the next 48 hours.
- Avoid exposing treated skin to direct sunlight for 2 hours after treatment. It is recommended to perform the treatment at night and, if you go out, apply sunscreen.
- Do not swim in pools for 48 hours after treatment.

COURSE OF TREATMENT

- **First month (once every 3–5 days):**

You will not notice any immediate changes in the appearance of your hair or skin. Treatment requires continued use.

- **Second month (once a week):**

The hair will become finer and softer. You will see the black hair follicles falling out during your shower.

- **Third month (once every 15 days):**

Most of the new hair will stop regenerating strongly and you can gently remove it. Only a little hair will remain in the area.

- **Fourth month (once every 30 days):**

The hair will disappear. You will have permanently removed hair from the treated area and achieved perfect skin.

By following this treatment guide, the device will help you permanently remove hair from specific areas. Some parts won't require any additional maintenance; others may need touch-ups every 6, 9, or 12 months.

REPLACING THE LAMP

1. Press and hold the power button for 2 seconds to turn off the device, then unplug the power cord.
2. Hold the top of the handle with one hand and the lamp with the other, pull it out and insert the corresponding one according to the function (hair removal - yellow lens, skin rejuvenation / acne treatment - red lens).

Lamp lifespan


Hair removal lamp (yellow lens): 300,000 pulses.

Skin Rejuvenation Lamp (Red Lens): 100,000 pulses.

SKIN REJUVENATION FUNCTION

1. Replace the skin rejuvenation lamp. The device will automatically test and identify the lamp type after power on.
2. Press the level adjustment button until it reaches level 1, direct the lamp to the area requiring skin rejuvenation to irradiate it.
3. It is recommended to irradiate each area once and perform the treatment once or twice a week. You will be able to observe a noticeable effect after maintaining the treatment for 2–3 months.

AREAS SUITABLE FOR TREATMENT

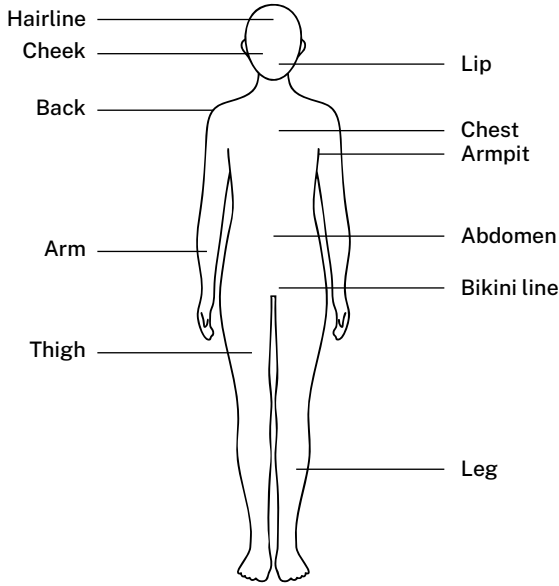ADVERSE REACTIONS AND RELATED COUNTERMEASURES

All cosmetic methods carry some risks, and this product is no exception. Possible adverse reactions and corresponding measures are:

- **Pain:** Since each person has different skin, the sensation of pain varies. It is rare for users to feel uncomfortable. A mild burning sensation lasting more than a minute may occur, but most adults and teens can tolerate it.
- **Mild skin redness or burning sensation:** It is a normal reaction. In most cases, this discomfort will disappear on its own after a short time. If it persists for more than 24 hours, please contact customer service.
- **Smell of burnt hair:** Indicates that you have not shaved all the hair before using the appliance. It's normal.

TECHNICAL SPECIFICATIONS

PARAMETER	WORTH	PARAMETER	WORTH
Power	36W	Cable length	1.8 m
Strain	Input:100-240V~ 50/60Hz Exit: 24V 1.5A	External environment	Temperature: -20~60°C Relative humidity: 20-90% non-condensing humidity
Weight	0.32 kg	Fluency	4.7 J
Irradiation area	4 cm ²	Dimensions	182x85x151 mm

TROUBLESHOOTING

PROBLEM	CAUSE	SOLUTION
No power	The power line is not properly connected	Connect the power line. Once connected correctly, the level indicator will light up for 2 seconds.
The device cannot be turned on	<ol style="list-style-type: none"> 1. The power button is not pressed correctly. 2. The power button is not held down for 2 seconds. 	When turning on, apply a little force and hold the power button for at least 2 seconds.
The level indicator flashes after switching on	The lamp is not installed correctly	Reinstalling the lamp: Turn off the appliance and unplug the power cord before replacing the lamp.
No light is emitted when the light trigger button is pressed	<ol style="list-style-type: none"> 1. The lamp is not in contact with the skin. 2. The lamp life has expired. 3. The charging process ended before the readiness indicator lit up. 	<ol style="list-style-type: none"> 1. Make the lamp touch the skin. 2. Replace the lamp. 3. Press the light trigger button after the readiness indicator is on.
Automatic shutdown	The device automatically switches off after 5 minutes if not in use.	Press and hold the power button for 2 seconds to restart.

In compliance with Directives: 2012/19/EU and 2015/863/EU on the restriction of the use of dangerous substances in electric and electronic equipment as well as their waste disposal. The symbol with the crossed dustbin shown on the package indicates that the product at the end of its service life shall be collected as separate waste. Therefore, any products that have reached the end of their useful life must be given to waste disposal centres specialising in separate collection of waste electrical and electronic equipment, or given back to the retailer at the time of purchasing new similar equipment, on a one for one basis. The adequate separate collection for the subsequent start-up of the equipment sent to be recycled, treated and disposed of in an environmentally compatible way contributes to preventing possible negative effects on the environment and health and optimises the recycling and reuse of components making up the apparatus. Abusive disposal of the product by the user involves application of the administrative sanctions according to the laws.

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