
CREATE

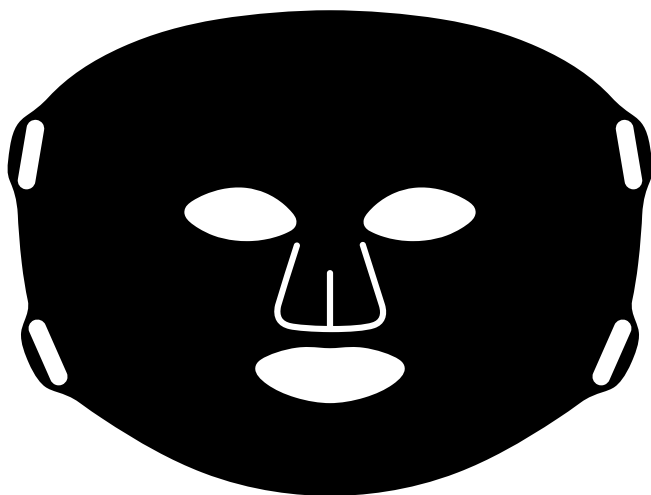Skin Care Led Mask

User manual | Manual de instrucciones

[illegible]

[illegible]

INDEX

ENGLISH

Safety warnings	6
Parts list	7
Command	7
Product introduction	8
Colors of light	8
How to use	9
Charging method	10
Cleaning and maintenance	10
Technical specifications	11

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ENGLISH

Thank you for choosing our product. Before using this device and to ensure its best use, please read the instructions carefully.


This manual contains safety information for the correct use of the product: read it carefully and retain it to avoid damage and hazards caused by improper use. The product is for personal use and should not be used by multiple people.

SAFETY WARNINGS

Before using any electrical appliance, carefully read the basic safety instructions.


- Do not use the device if you are particularly sensitive to light, unless a sensitivity test has been performed.
- Do not use the device if you are pregnant, as the risks are unknown and have not been established.
- Do not use the device if you suffer from epilepsy, as some lights or frequencies may trigger a seizure.
- The device is not waterproof, do not allow the entry of water or other liquids as this may impair its operation.
- Do not use the device if you are taking medication that may cause skin sensitivity to light.
- Do not use the device while sleeping, walking, or driving.
- Use the device according to the correct instructions and avoid excessive use.
- Discontinue use if discomfort occurs or if skin redness or discoloration lasts more than 24 hours. Consult your dermatologist.
- Do not touch the LED bulbs while the device is turned on.
- Do not use the device if it is visibly damaged.
- Do not subject the device to impacts or shocks to avoid deformation or damage.

PARTS LIST

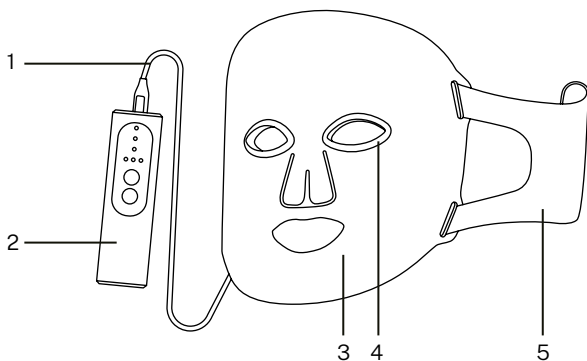- | | |
|---------------------------|------------------|
| 1. USB cable | 4. Eye masks |
| 2. Command | 5. Elastic bands |
| 3. Phototherapy face mask | |

COMMAND

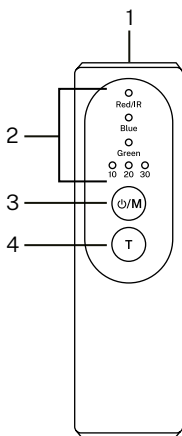- | |
|--------------------------|
| 1. USB-C port |
| 2. Mode indicators |
| 3. Power/light button |
| 4. Time selection button |

PRODUCT INTRODUCTION

This device uses the principle of phototherapy to activate skin cells and improve their condition by releasing specially processed light at specific wavelengths. This is a non-invasive, chemical-free, painless, and safe skin care method.

The light, once absorbed by the skin, is transformed into intracellular energy, triggering a photochemical and enzymatic reaction. During this process, a large amount of glycogen and proteins are produced, which promote cellular metabolism and collagen regeneration, thus achieving whitening, rejuvenation, acne treatment, spot removal, skin firming, and wrinkle reduction effects.


COLORS OF LIGHT

The device uses lights in 8 different wavelengths, with the following main functions:


1. **Red light (630 nm $\pm$ 10 nm; indicator light in red color):** Called bioactive light, it is characterized by its high purity, intensity and uniform density. It is effective in skin care and health, as it increases cellular vitality, accelerates blood circulation, boosts metabolism, and promotes the regeneration of collagen and fibrous tissue. Thanks to this, it also improves the elasticity and dullness of the skin, being indicated for anti-aging, antioxidant, and repair effects.
2. **Infrared light (850 nm $\pm$ 10 nm, invisible; indicator light in dark red):** It promotes collagen production, accelerates skin metabolism, brightens skin tone, reduces pores, and provides a more youthful appearance.
3. **Blue light (450 nm $\pm$ 10 nm; indicator light in blue color):** It coincides with the absorption peak of porphyrin, a metabolite of *Propionibacterium acnes*. During the de-excitation process following stimulation, active oxygen is generated, creating a highly oxidative environment capable of eliminating bacteria, improving and treating acne.
4. **Green light (520 nm $\pm$ 10 nm; indicator light in green):** It helps soothe the skin, relax the nerves, restore balance to sensitive skin, stabilize and relax the nervous system, and promote the elimination of toxins.
5. **Yellow light (580 nm $\pm$ 10 nm; combination of red and green indicator lights):** It improves oxygenation of skin cells, promotes microcirculation, breaks down and attenuates pigmentation, stimulates cell growth, softens the skin, brightens its tone, repairs the dermis, and reduces wrinkles.

6. **Purple light ($420\text{ nm} \pm 10\text{ nm}$; combination of red and blue indicator lights):** It has excellent healing effects in treating and removing acne marks.
7. **Cyan light ($490\text{ nm} \pm 10\text{ nm}$; combination of blue and green indicator lights):** Increases cellular energy and promotes metabolism.
8. **White light ($400\text{--}700\text{ nm} \pm 10\text{ nm}$; combination of red, blue, and green indicator lights):** It penetrates the deep layers of the skin, accelerates metabolism, breaks down pigmentation, and improves lines and sagging skin.



HOW TO USE

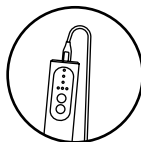1. Wash your face and let it air dry naturally or pat it dry gently with a soft towel.
2. Take off the mask. Insert the elastic bands into the corresponding slots and connect the cable to the controller.

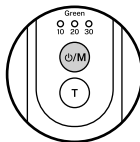3. Place the mask over your face. Hold both sides of the elastic bands and adjust the Velcro at the back of your head. Adjust the straps to ensure a secure fit and prevent light leakage through the silicone eye protector.
4. Press and hold the power button for 3 seconds to turn on the mask. The default mode is red light mode.

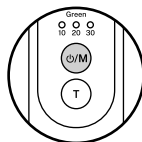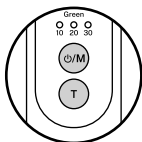5. Press the power button to select the mode and the time selection button to choose the duration. There are 8 modes: red, infrared, blue, green, yellow, purple, cyan and white light. During use, you can switch modes and select a 10, 20, or 30-minute timer.
6. The device will automatically turn off when the time is up. You can also press and hold the power button for 3 seconds to manually turn it off. Note: It is recommended to use the device 3 to 5 times a week.

NOTE: Since infrared light is a type of invisible light, if you select this mode you will not see an obvious bright light, just a slight reddish tint to the LEDs when viewed closely. This is a normal phenomenon.

CHARGING METHOD

1. Connect the controller with an approved adapter (not included) using the Type C charging cable.
2. Plug it into the power outlet. The indicators will light up and the battery symbol will scroll during charging.
3. When the battery symbol remains solid, it means the charge is complete. Then disconnect the adapter.

CLEANING AND MAINTENANCE

- Clean the device with a dry cloth after each use. Do not rinse under water.
- Store in a cool, shaded place, avoiding direct sunlight or high temperatures.
- Keep the device out of reach of children to avoid risks.
- Do not attempt to open or repair the device yourself. If any failure occurs or the product stops working, please contact customer service.
- Store the device in the storage bag after use.

TECHNICAL SPECIFICATIONS

PARAMETER	WORTH
Power	4.2 W
Strain	5 V
Intensity	1 A
Battery	2600 mAh
Cable type	USB-C
Loading time	3/4 h
Operating time	120 minutes approx (maximum power), 240 minutes approx (minimum power)
Dimensions	Mask: 284x217x20 mm Control: 37x27x129 mm
Weight	Mask: 193 g Control: 89 g
Materials	Silicone

In compliance with Directives: 2012/19/EU and 2015/863/EU on the restriction of the use of dangerous substances in electric and electronic equipment as well as their waste disposal. The symbol with the crossed dustbin shown on the package indicates that the product at the end of its service life shall be collected as separate waste. Therefore, any products that have reached the end of their useful life must be given to waste disposal centres specialising in separate collection of waste electrical and electronic equipment, or given back to the retailer at the time of purchasing new similar equipment, on a one for one basis. The adequate separate collection for the subsequent start-up of the equipment sent to be recycled, treated and disposed of in an environmentally compatible way contributes to preventing possible negative effects on the environment and health and optimises the recycling and reuse of components making up the apparatus. Abusive disposal of the product by the user involves application of the administrative sanctions according to the laws.

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